

Product Datasheet


Anti-DARPP-32




Overview

Catalog #	MADARPP32
Host Species	Rabbit Polyclonal
Format	Antigen Affinity Purified from Pooled Serum
Applications	WB 1:1000 IHC 1:250
Species Tested	Mouse, Rat
Expected Reactivity	Bovine, Canine, Chicken, Human, Non-Human Primate, <i>Xenopus</i> , Zebrafish
Immunogen	Synthetic peptide corresponding to amino acid residues from the N-terminal region of rat DARPP-32, conjugated to keyhole limpet hemocyanin (KLH).
Molecular Weight	32 kDa
Cite this Antibody	PhosphoSolutions Cat# 374-DARPP, RRID:AB_2492067

Images

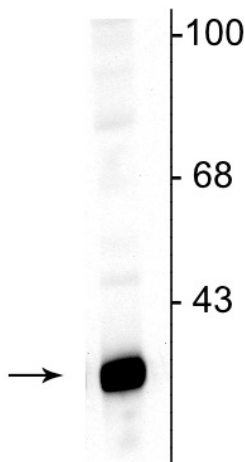

Western blot of rat hippocampal lysate showing specific immunolabeling of the ~32 kDa DARPP protein.

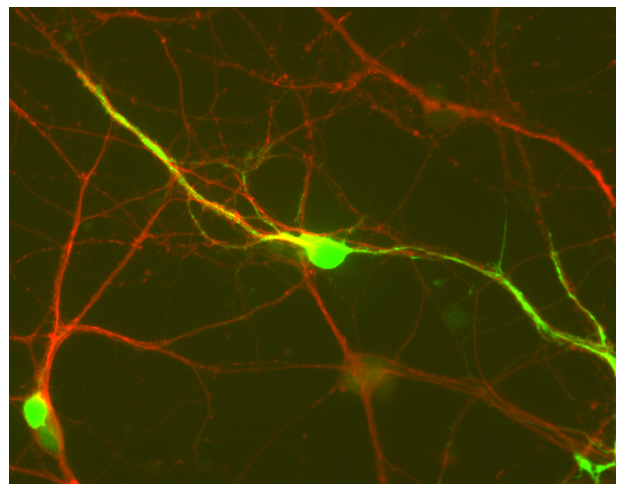

Immunochemical staining of cultured CD1 mouse striatal cells showing specific labeling of DARPP-32 (Cat. # 374-DARPP, green, 1:250) and Tuj1 (red). Cells and photo courtesy of QBM Cell Science.

Details

Target Description	DARPP-32 is a dopamine (DA) and cAMP-regulated ~32 kDa phosphoprotein that is associated with dopaminergic neurons (Fienberg et al., 1998). The protein inhibits protein phosphatase I when it is phosphorylated on Thr-34. In contrast, when DARPP-32 is phosphorylated on Thr-75 the protein acts as an inhibitor of PKA (Bibb et al., 1999). Phosphorylation of DARPP-32 is thought to play a critical role in the regulation of dopaminergic neurotransmission. In addition, the activity of DARPP-32 is also thought to play important roles in the actions of alcohol, caffeine and Prozac® (Maldve et al., 2002; Lindskog et al., 2002; Svenningsson et al., 2002).
Specificity	Specific for endogenous levels of the ~32 kDa DARPP-32 protein.
Production/Purification	Prepared from pooled rabbit serum via chromatography on an affinity column prepared with the N-terminal peptide used as antigen.
Quality Control	Western blots performed on each lot.
Buffer	10 mM HEPES (pH 7.5), 150 mM NaCl, 100 µg per ml BSA and 50% glycerol.
Storage	Storage at -20°C is recommended, as aliquots may be taken without freeze/thawing due to presence of 50% glycerol.
Stability	After date of receipt, stable for at least 1 year at -20°C.

Significant Citations

Bonito-Oliva, A., DuPont, C., Madjid, N., Ögren, S. O., & Fisone, G. (2015). Involvement of the striatal medium spiny neurons of the direct pathway in the motor stimulant effects of phencyclidine. *The international journal of neuropsychopharmacology (CINP)*.

Michelle Niculescu, Shane A. Perrine, Jonathan S. Miller, Michelle E. Ehrlich, and Ellen M. Unterwald (2008) Trk: A Neuromodulator of Age-Specific Behavioral and Neurochemical Responses to Cocaine in Mice. *J. Neurosci.*, Jan 2008; 28: 1198 - 1207.

Our Guarantee

As an original manufacturer, we are dedicated to creating quality and reproducible antibodies that further your research. We provide personalized customer support from the scientists that made the antibody and offer a free replacement or 100% refund if we cannot resolve an issue. Order today and experience **Antibodies that Work™**.

Note: For research use only. Not intended for therapeutic or diagnostic use. Use of all products is subject to our terms and conditions, viewable on our website.

